CeQur Simplicity™

Mealtime Insulin Patch

1U per click

IMPORTANT REMINDERS

- Please read all instructions before using your CeQur Simplicity™ Patch (or 'Patch').
- Use the Patch only with U-100 rapid-acting (mealtime) insulin: Humalog® or NovoLog®.
- Use your Patch for bolus insulin needs at mealtimes and snack times, or for correcting high blood glucose.
- Do not fill your Patch with long-acting insulin or any other insulin or medication.
- Do not directly expose your Patch to oils, lotions, or sunscreen.
- Do not wear your Patch for longer than 7 days.

Refer to User Guide for full instructions by scanning the QR code.


*Visit **myceqursimplicity.com** for more information.*

CeQur Simplicity™

Mealtime Insulin Patch

1U per click

Each click delivers **1 unit** of insulin

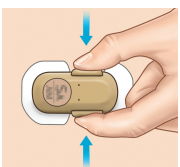

IMPORTANT

- The first click of buttons for a new patch delivers 0.3 units short of 1 unit.
- If buttons lock during a dose, the last click of the buttons may deliver less than 1 unit.
- If uncertain about your delivered dose, stop and check your blood glucose in 1 to 2 hours.

©2026 CeQur Corporation

AW-MTL-00877 Ver. 1